


## Therma-Liner™ Bag Preparation

- Use one Liner Bag per single patient only (12-hours)
- Record information on the Liner Bag (patient name, care provider name, date and time, other shift information)
- Open outer water compartment (blue zipper) and fill with room temperature water up to the fill mark
- To close, press firmly along the entire blue zipper for a tight seal (Liner Bag prep required only once per 12-hour use)
- Carefully insert feeding container into the inner compartment of the Liner Bag. Make sure to gently push the container to the bottom of the Liner Bag (Note: Use care when removing and re-inserting milk containers from the Liner Bag to prevent any potential damage)
- Carefully place the Liner Bag with feeding container into the well of the warmer. Make sure to gently push to the bottom of well.



## To Warm/ Thaw Feeding


- Turn the device on by flipping the rear power switch ON (standby mode). Then press the on/off button to begin.
- Scroll up or down to choose desired mode from the below three options and press button corresponding to the green check mark to select:

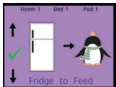

Fridge To Feed



Frozen To Feed

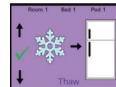

Thaw

- Choose the closest match for container / volume size  
*At this point the Penguin will flash a message on LCD screen to insert the Therma-Liner Bag. If Liner Bag has already been inserted per Therma-Liner Bag preparation instructions above proceed to the next steps. If Liner Bag has not been inserted yet, follow Therma-Liner Bag preparation instructions above to continue)*

HOLD time is provided for warm modes only. For Thaw mode, remove feeding immediately

- Press ENTER to begin Warming
- The screen turns GREEN once the feeding is ready. You will also hear a beeping sound
- After a warming cycle is complete, there is 30 minutes of Hold time, then the screen turns red. The screen turns RED once it's past the 30 minutes of HOLD time after warming cycle is complete.
- Remove the recyclable Therma-Liner Bag and press ON/OFF button to turn off

